USER MANUAL

ENSURE YOUR BEST RESULTS

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3) Energizing Hum Indicates The Device Is On (Not Shown)

Un bourdonnement énergisant indique que l'appareil est allumé (non illustré)
El zumbido energizante indica que el dispositivo está encendido (no se muestra)
Brummen zeigt an, dass das Gerät eingeschaltet ist (nicht abgebildet)
Ronzio stimolante che indica che il dispositivo è acceso (non mostrato)
윙윙거리는 소리는 기기가 켜져 있음을 나타냅니다(표시되지 않음)
电蜂鸣声表明仪器已开启 (未显示)

**4) White LED Power/
Treatment Modes**

Voyants DEL blancs d'alimentation/intensité
Luces led de encendido/intensidad blancas
Weiße LED-Leuchten = Power und Intensität
Luci di accensione LED bianche/Luci di intensità
백색 LED 전원/강도 표시등
白色LED功率/强度灯

**5) NuFACE FIX® Line Smoothing Serum
(Microcurrent Activator) 0.5 fl oz / 15 mL**

Sérum NuFACE FIX® (activateur de microcourant) 0,5 fl oz/15 mL
Sérum lisador de líneas NuFACE FIX® (activador de microcorrientes) 0,5 fl. oz./15 mL
NuFACE FIX® Glättendes Serum (Mikrostromaktivator) 0,5 fl oz / 15 ml
Siero levigante NuFACE FIX® Line (Attivatore di microcorrente) 0,5 fl oz / 15 mL
NuFACE FIX® 라인 스무딩 세럼 (마이크로커런트 액티베이터) 0.5 fl oz / 15 mL
NuFACE FIX®平滑細紋精華 (微电流激活器) 0.5 fl oz / 15 mL

**6) USB Charging Cable
(Not Shown)**

Câble de chargement USB (non illustré)
Cable de carga USB (no se muestra)
USB-Ladekabel (nicht abgebildet)
Cavo ricarica USB (non mostrato)
USB 충전 케이블 (표시되지 않음)
USB充电线 (未显示)

**7) USB to AC Wall Adaptor
(Sold Separately/Not Shown)**

Adaptateur mural USB vers CA (vendu séparément/non illustré)
Adaptador de USB a toma de corriente CA (se vende por separado/no se muestra)
USB-zu-AC-Wandadapter (separat erhältlich/nicht abgebildet)
Adattatore a parete da USB a CA (Venduto separatamente/Non mostrato)
USB-AC 벽면 어댑터 (별도 판매/표시되지 않음)
USB转AC墙式适配器 (单独销售/未显示)

**N U F A C E®
FIX MicroWand**

The FIX MicroWand device produces low levels of microcurrent which is emitted through the two spherical probes as they both touch the skin. Prior to using the FIX MicroWand device, the selected area of the skin to be treated is first prepped with a NuFACE® Microcurrent Activator to moisturize the skin and ensure optimal microcurrent conductivity. The two spherical probes are then gently glided over the skin to deliver the electrical impulses per the selected treatment setting. Microcurrent is an aesthetic modality providing electrical current in millionths of an ampere and has the ability to increase the firmness of the skin and muscles.

INTENDED USE

The NuFACE FIX MicroWand Device is intended for facial stimulation and is indicated for over-the-counter cosmetic use. Any consequence resulting from misuse, application of the device to other parts of the body or face that are not approved for use, connection to improper voltage sources, dirty, spares, expired or non-NuFACE conductive alternative, or any other improper application is neither the responsibility of NuFACE nor its affiliates.

DO NOT use on infants and children under the age of 18. Keep this device out of the reach of children and pets.

IMPORTANT INFORMATION

- Use the FIX MicroWand Device only as described within this User Manual.
- The risks and dangers of using the device in any way that does not follow this User Manual are unknown.
- Your NuFACE FIX MicroWand Device has a unique serial number located on the back of the device. Record your serial number and register your device online at www.NuFACE.com.
- The NuFACE device is designed for cosmetic use only.
- The device should be operated, transported, and stored in the manner specified in this manual.
- Failure to use and maintain the device and corresponding attachments in accordance with the instructions outlined in this manual will void the warranty and may result in product damage or physical injury.
- The long-term effects of transcutaneous electrical stimulators are unknown.
- Apply stimulation only to normal, intact, clean, healthy skin.
- Avoid applying excessive pressure when using the device.
- Skin irritation or hypersensitivity may occur due to the electrical stimulation and / or use of the NuFACE® Microcurrent Activator.
- A microcurrent sensation may occur while using the device. Using a lower intensity and/or additional activator may be needed to reduce the sensation.
- Skin may become dry during treatment in which Microcurrent Activator may need to be reapplied.
- A perception of flashing lights may occur if device is used near or around the eye region due to stimulation of retina and/or optic nerve from the electrical current. The perception of flashing lights is not harmful.
- **DO NOT** use the device if you see any signs of damage or leaks from the device.
- **DO NOT** remove the battery of the device.
- **DO NOT** store or leave the device under the sun or on a hot surface.
- Follow all local regulations and waste management guidelines for safely disposing devices containing batteries.
- The device has been designed to provide freedom from unacceptable risk when used under normal conditions and reasonably foreseeable hazards. In the event of complete loss of treatment function, the device does not pose a risk to the user.

CONTRAINDICATIONS FOR USE

Consult your healthcare provider if you have any pre-existing or potential health conditions that may be impacted by use of the device including those listed below.

- **DO NOT** use if you have an electronic implanted device such as a defibrillator, neurostimulator, pacemaker, or ECG Monitor. If you use a body-worn medical device, such as an insulin pump, consult your healthcare provider.
- **DO NOT** use this device if you have impaired cognition.

CONTRAINDICATIONS FOR USE (CONTINUED)

- **DO NOT** use this device if you have an infection, tuberculosis, or swelling of bone tissue.
- **DO NOT** use if you are pregnant.
- **DO NOT** use if you have any active cancer or over any suspicious or cancerous lesion.
- **DO NOT** use this device if you have a history of epilepsy or seizures.
- **DO NOT** use this device if you have a history of cardiac failure.
- **DO NOT** use in areas of impaired skin sensation
- **DO NOT** use if you have hemorrhagic conditions or bleeding disorders
- **DO NOT** use this device on any area of the body other than instructed in this manual including across the chest, the temples, upper back, and groin area.
- **DO NOT** use this device directly on the anterior neck, over the carotid artery/sinus.
- Consult your healthcare provider if any of the contraindications above apply to you.

PRECAUTIONS

- If you are under the care of a physician, consult your physician before using this device.
- Apply stimulation only to normal, intact, clean, healthy skin.
- **DO NOT** use Microcurrent Activator or FIX MicroWand directly on lips or eyes.
- Use caution when using this device near metal implants, locally over areas of swollen, infected, or inflamed areas or skin eruptions, e.g., phlebitis, thrombophlebitis, broken capillaries, varicose veins, etc.
- Use caution when using this device if you have deep-vein thrombosis, impaired circulation, chronic wounds, or skin disease (e.g. Eczema)
- Use caution to avoid applying electrical stimulation near the thyroid gland.
- **DO NOT** operate FIX MicroWand Device in the presence of electronic monitoring equipment (e.g., heart monitors, heart tracing or in close proximity (eg. 1 m) to shortwave, microwave or electromagnetic therapy equipment as it may cause the medical equipment to malfunction.
- Use of accessories and cables other than those sold by NuFACE® may result in increased electronic emissions, decreased electronic immunity of the system, damage to the device, or injury.
- Only use the external charging cable provided with the FIX MicroWand Device. Optional SVAC 1A USB to AC wall adapters can purchased directly from NuFACE or authorized NuFACE sellers. The optional power adapter is marked NuFACE®. Always position the power adapter so it is accessible in case it is needed to be disconnected from the wall outlet.

⚠ WARNINGS

DO NOT apply stimulation or use the device:

- Directly on lips or eyes.
- While bathing or showering.
- At extreme temperatures.
- If you suspect the device is damaged or not functioning properly.
- While driving or operating machinery.
- **DO NOT** disassemble the device as this may cause damage, malfunction, electrical shock, fire, or injury. There are no user-serviceable parts inside.

Note: Disassembly voids the warranty

Discontinue use of the device and contact your healthcare provider if you experience any of the following:

- Persistent pain or discomfort.
- Discomfort or any difficulty breathing or swallowing.
- Negative or harmful reaction(s) or injury.
- Symptoms of hypertension (dizziness or lightheadedness).
- Constant headaches and/or painful sensations.
- Skin irritation or burning.
- Unpleasant sensations on skin during use.
- Consistent flashing lights without the use of your device.
- If you experience any negative reactions while using the device and/or the microcurrent activator.

BEFORE USING YOUR NEW DEVICE

Follow these important steps before using your new device:

- 1) Remove from the box your FIX MicroWand device, USB charging cable.
- 2) Connect your FIX MicroWand Device to a USB outlet to begin charging, safely away from any source of water, on a counter or tabletop. The lights will increasingly illuminate, indicating that your device is charging. All of the lights will stay illuminated when device is fully charged.


Note: The FIX MicroWand Device will not turn on when it is charging.

- 3) **IMPORTANT: Charge your FIX MicroWand Device for 4 hours before the first use.** The white LED light flashes to indicate that the beauty device is charging. When charging is complete, the white lights will light up continuously.
- 4) **IMPORTANT: Test your skin to make sure you are not sensitive to the NuFACE Microcurrent Activator.** To do so, rub a dime-sized amount of Microcurrent Activator on your forearm and leave for 15-30 minutes before washing it off. If an allergic reaction (rash, redness, etc) develops, immediately contact NuFACE for recommendations on alternative NuFACE formulas. You should not initiate treatment with the FIX MicroWand Device until your reaction to the Microcurrent Activator has been resolved.
- 5) When your FIX MicroWand Device is ON it will emit an energizing hum and the white LEDs will illuminate.
- 6) **IMPORTANT: Always ensure that the two (2) spheres remain in direct contact with the skin.** Avoid pulling or stretching the skin during treatment.
- 7) The FIX MicroWand Device energizing hum will emit 3 rapid pulses during treatment when the battery needs recharging. Note: The FIX MicroWand Device will NOT contract the muscles. The current is soft, soothing, and mild when used correctly.

TREATMENT AREAS


Help address fine lines and wrinkles around these treatment areas: Forehead, Eye Area, and Mouth Area (Nasolabial Folds).

IMPORTANT: Avoid using the device directly on your lips or eyes.

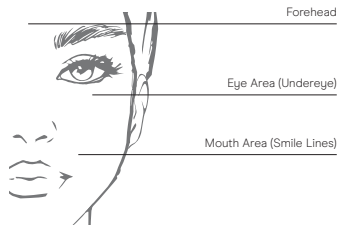TREATMENT TECHNIQUES

Ensuring that both spheres are in contact with the skin, apply light to medium pressure while using the FIX MicroWand Device feathering technique.

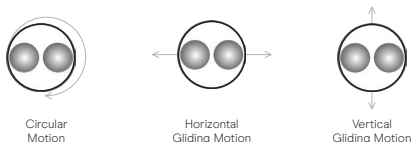DEVICE INTENSITY INDICATOR

LED LIGHT INDICATOR	TREATMENT MODE
1 LED	●○○○ Eye Area (Undereye)
2 LED	●●○○ Mouth Area (Smile Lines)
3 LED	●●●○ Fine Lines (Forehead)

HOW TO USE YOUR DEVICE

The instructions below are designed to provide the most effective results when followed correctly.

1) CLEANSE

- Cleanser your face with an oil-free cleanser. Men should shave before treatment — facial hair can interfere with the conductivity of the spheres.

Note: Any cream, lotion, self-tanner, and/or oil-based products on the skin can block the microcurrent.

2) ACTIVE

IMPORTANT: Always apply a NuFACE Microcurrent Activator when using your device.

- Apply a mask-like layer of the NuFACE Microcurrent Activator section-by-section as you treat.
- Reapply activator if skin becomes dry during treatment in order to continue safely using device.

3) PREPARE YOUR DEVICE

- Unplug the USB Cable from your FIX MicroWand Device and turn it ON by pressing the Power Button; the white LEDs will illuminate, and the energizing hum will begin when the device is ON.

4) BEGIN TREATMENT

- I. Pick an area on your face to begin your treatment (See **TREATMENT AREAS**)
- II. The FIX MicroWand device has three treatment modes. The treatment mode can be changed with a quick press of the Power Button.
- III. Select the proper LED Light Indicator/Treatment Mode (See **DEVICE INTENSITY INDICATOR**)
- IV. Place the device directly on your skin and keep both spheres in contact during treatment to keep the microcurrent flowing.
- V. Apply light to medium pressure to begin treatment.
- VI. During treatment, apply any of the following techniques: Circular, Horizontal, Vertical motions (See **TREATMENT TECHNIQUES**).

You can treat one or more treatment areas for the full 3-minute cycle.

IMPORTANT: If you feel the electrical current and / or experience a tingling sensation on your skin, do the following:

- **Add more NuFACE Microcurrent Activator to the area you are treating.** The device will NOT contract the muscles when used as instructed.

5) FINISH TREATMENT

- Your FIX MicroWand Device will automatically turn OFF after a 3-minute treatment.
- The device can also be turned OFF by holding the power button while the device is running until the indicator LEDs turn off, and the device stops humming.
- After treatment is complete and device is turned off, massage remaining NuFACE Microcurrent Activator into skin followed by eye cream, moisturizer, and SPF.

Go to the section '**CARING FOR YOUR DEVICE**' for instructions on cleaning and maintaining your device.

TREATMENT FREQUENCY

- The FIX MicroWand Device may be used up to 2 times a day per treatment area.
- 5 days a week.
- Can be used as the finishing step to your NuFACE TRINITY+ or MINI+ treatment.

For more advanced treatment options and holds, download the **NuFACE Smart App**, or visit our website at www.NuFACE.com.

RESULTS

Results may vary based on skin type and skin condition. Optimal results are achieved with regular use.

CARING FOR YOUR FIX MicroWand DEVICE

Storage:

The device is intended to be stored and used at room temperature. Use caution prior to using the device if it has been stored or not used for an extended period of time. Inspect the device and do not use it if there are any signs of internal battery leaks.

Cleaning:

• Slightly dampen a soft, lint-free cloth with water or rubbing alcohol to clean your device after every treatment. Visually inspect the cleaned surfaces for contamination. If dirt or contamination is still present, repeat the cleaning steps described above. If contamination cannot be removed contact your call center for additional instructions.

- **DO NOT submerge your device in water at any time.**

Charging:

- Before first use, charge your FIX MicroWand Device for 4 hours.
- The device will NOT turn on when it is charging.
- For continued optimal battery life, it is recommended to charge your device regularly.
- After being fully charged, the device shall operate for up to 40 full 3-minute treatment cycles.
- Keep the product out of the reach of children when charging.
- Keep the product away from high temperature, open flame and water sources when charging to prevent short circuiting or damaging the battery.
- Plug the charger into a wall outlet away from water, and within a location that can be accessed and disconnected easily.

Battery Care:

For ideal battery life, we recommend charging your device when the low battery indicator (1 LED light only) is displayed on your device and the low battery indicator begins.

Check Battery:

To check your device battery level, **hold the power button down on your device for one (1) second**. The battery level will be illuminated through the three (3) LEDs:

LED LIGHT INDICATOR	BATTERY POWER
No LED Lights	○○○ No Battery Power A full charge cycle of four (4) hours is recommended before use.
1 LED Flashing	⦿○○ Critical Battery. Not enough charge to start a treatment.
1 LED	●○○ Low Battery
2 LEDs	●●○ Medium Battery
3 LEDs	●●● Full Battery

CARING FOR YOUR FIX MicroWand DEVICE (CONTINUED)

Note: The rechargeable battery of the FIX MicroWand Device is not removable by the user. Use caution prior to using the device if it has been stored or not used for an extended period of time. Inspect the device and do not use it if there are any signs of internal battery leaks.

Disposal:

IMPORTANT: If disposing of the device, do not put into a fire. Dispose of the device and its batteries per local regulations.

TROUBLESHOOTING

Device Not Turning On:

Check to see if your device is fully charged.

- If the device does not turn ON, the battery may be drained.
- A full charge cycle of four (4) hours may be required for the device.
- For additional information go to the "Battery" section in this user guide.

Resetting Your Device:

- To reset your device, plug your device into the power supply.
- Hold down the power button for **10 seconds**.
- The 1st and 3rd LEDs will blink three (3) times.
- Once the LEDs stop blinking you may start using your device.

DO NOT attempt to troubleshoot by disassembling the product as this will void the warranty and may cause injury.

CUSTOMER SUPPORT

If you are still experiencing issues with your device, contact NuFACE Client Services at info@nuface.com or call **888-360-9171** (US Only).

WARRANTY

NuFACE® warrants that all NuFACE® devices and attachments are free from manufacture defects in materials and workmanship for a period of two (2) years from the date of your original purchase. Our obligation under this Warranty is limited to replacement of the product which will not extend the period of your Warranty. This Warranty does not cover damage covered by abuse, improper care, application of harmful chemicals, product that may have been altered in any way, loss of the device and product purchased from an unauthorized seller.

NO OTHER WARRANTY EXPRESSED OR IMPLIED, WHETHER OF FITNESS FOR A PARTICULAR PURPOSE OR OF MERCHANTABILITY OR OF ANY OTHER KIND, SHALL EXIST IN RESPECT TO SUCH PRODUCT.

For Devices Purchased in the US: For more information about your warranty or to initiate a warranty claim and receive a return authorization number, contact our Client Services Team via our website at www.NuFACE.com or by calling **888-360-9171**. A receipt must be shown from an authorized seller as proof of purchase.

For Devices Purchased Outside of the US, International Returns, and International Warranties: For all product purchased outside of the United States, warranty and returns are per your place of purchase—spa, retailer, or distributor. Please consult your receipt or place of purchase for their return and warranty policy.

TECHNICAL INFORMATION

The NuFACE® FIX MicroWand complies with the requirements of IEC 60601-1-2:2014 (EMC Collateral Standard) including the E-field susceptibility requirements at a 10 volts per meter, at frequencies from 80 MHz to 2.7 GHz. However, even at this level of device immunity, certain transmitting devices (cellular phones, two-way radios, cordless phones, paging transmitters, RFID devices, etc.) emit radio frequencies that could interrupt FIX MicroWand operation if operated in a range too close to the device. Users should be aware of possible radio frequency interference if portable devices are operated in close proximity to the FIX MicroWand. It is recommended that transmitting devices such as cellular phones, two-way radios, cordless phones, etc. not be located within 10 feet of the FIX MicroWand during treatment. The NuFACE® FIX MicroWand complies with the requirements of IEC 60601-1-2:2014 for Radiated and Conducted Emissions in Group 1 Class B.

The NuFACE® FIX MicroWand complies with the requirements of IEC 60601-1-2:2014 for Harmonic distortion in Harmonic classification A.

ELECTROMAGNETIC COMPATIBILITY INFORMATION

- The FIX MicroWand device is intended and suitable for use in a Home Healthcare Environment.
- The FIX MicroWand device does not incorporate any form of radio frequency or wireless technologies and was tested to not cause harmful interference to other devices.
- The FIX MicroWand functions and performance are not expected to be degraded when operated within a Home Healthcare Environment.
- When cleaning the device, follow the care and maintenance instructions provided in this User Manual.
- When charging the device, follow the instructions provided in this User Manual.

The FIX MicroWand device was tested to be compliant with the applicable Electromagnetic Compatibility (EMC) tests per IEC 60601-1-2 (2020) and IEC 60601-4-2.

Emissions

TEST	TEST STANDARD	COMPLIANCE CLASSIFICATION
Radiated RF Emissions	CISPR 11	(Group 1, Class B) Home Healthcare Environment
Conducted RF Emissions Flashing	CISPR 11	(Group 1, Class B) Home Healthcare Environment
Harmonic Distortion	IEC 61000-3-2	Compliant
Voltage Fluctuations and Flicker	IEC 61000-3-3	Compliant

Enclosure Immunity

TEST	TEST STANDARD	COMPLIANCE CLASSIFICATION
Electrostatic Discharge (ESD)	IEC 61000-4-2	± 8kV contact ±2 kV, ±4 kV, ±8 kV, ±15kV air
Radiated RF EM Field	IEC 61000-4-3	10 V/m 80 MHz to 2.7 GHz 80% AM at 1 kHz
Proximity Fields from RF Wireless communications Equipment	IEC 61000-4-3	See RF Wireless Test Level table below
Rated power frequency magnetic fields	IEC 61000-4-8	30A/m 50 Hz or 60 Hz

ELECTROMAGNETIC COMPATIBILITY INFORMATION (CONTINUED)

Inport AC Power Immunity

TEST	TEST STANDARD	IMMUNITY TEST LEVEL
Electrical fast transient / bursts	IEC 61000-4-4	±2 kV 100 kHz repetition frequency
Surges Line-to-line	IEC 61000-4-5	±0.5 kV, ±1 kV
Surges Line-to-ground	IEC 61000-4-5	±0.5 kV, ±1 kV, ±2 kV
Conducted disturbances induced by RF fields	IEC 61000-4-6	3 V 015 – 80 MHz 6 V in ISM and amateur radio bands between 015 MHz and 80 MHz 80 % AM at 1 kHz
Voltage dips	IEC 61000-4-11	0 % UT; 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% UT; 1 cycle and 70% UT; 25/30 cycles Single Phase: at 0°
Voltage interruptions	IEC 61000-4-11	0% UT; 250/300 cycle

RF Wireless Test Levels

RF WIRELESS TEST LEVELS PER IEC 60601-1-2		
FREQUENCY RANGE (MHz)	SERVICE	IMMUNITY TEST LEVEL (V/M)
385	TETRA 400	27
450	GMRS 460, FRS 460	28
710	LTE Band 13, 17	9
745		9
780		9
810		GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5
870	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	28
930		28
1720		28
1845	28	28
1970		28
2450	Bluetooth, WLAN, 80211 b/g/n, RFID 2450, LTE Band 7	28

ELECTROMAGNETIC COMPATIBILITY INFORMATION (CONTINUED)

RF WIRELESS TEST LEVELS PER IEC 60601-1-2		
FREQUENCY RANGE (MHz)	SERVICE	IMMUNITY TEST LEVEL (V/M)
5240	WLAN 802.11 a/n	9
5500		9
5785		9

Conforms to EN IEC/UL: Device - 60601-1, 60601-1-6, 62366-1, 60601-1-11, 60601-2-10, 60601-1-2, 60529, 62133-2, 60335-1, 60335-2-23 Power Supply - 60950-1, 60068-2-70

SYMBOLS

One or more of the following markings may appear on the device, User Manual or packaging.

-  European Conformity
-  UK Conformity Assessed
-  Nemko Electrical Safety Certification Program
-  European Representative
-  Part Number
-  Operator's Manual; Operating Instructions
-  Refer to Instruction Manual/Booklet
-  Manufacturer
-  Complies with European Restriction of Hazardous Substances Directive 2002/95/EC
-  Product must be disposed of in accordance with the WEEE directive
-  Type BF applied part, which is equipment that provides a particular degree of protection against electric shock. Particularly regarding allowable leakage current and of the protective earth connection.
-  Class II double insulated electrical device
-  Protection against dust: Protection against splashed water
-  Detachable power supply part

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FIX MicroWand DEVICE GENERAL SPECIFICATIONS

- Device Measurements: 154mm/6.06"(L) x 27mm/1.06"(W) x 24mm/0.94"(D)
- Net Weight: 56 g / 2 oz
- System Kit Gross Measurements:
 - 197mm/7.75" (H) x 89mm/3.5" (W) x 59mm/2.34" (D)
- System Kit Gross Weight: 108 g / 5.8 oz
- Antenna Peak gain: 2.21 dBi
- Antenna type: Chip Antenna
- Antenna impedance: 50Ω nominal

Environmental Conditions

- Operating conditions: 10 to 40 °C, 30 to 75%rH*, 700 to 1060HPa
- Storage conditions: -10 to 55 °C, 10 to 90%rH*, 500 to 1060HPa
- Shipping condition: -10 to 55 °C, 10 to 90%rH*, 500 to 1060HPa

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- Battery Rating: 3.7V, 350mAh
- Max Output Current: 225µA
- Max Output Voltage: 28VDC
- Pulse Duration: 10ms – 60ms
- Pulse Frequencies: 8.3 – 50 Hz
- NuFACE AC Power Adapter DZ006THU050100H: 14030
- NuFACE USB Cable (5 V) to Barrel KC-DZ081A: 10671

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All Rights Reserved. / Tous Droits Réservés. / Todos Los Derechos Reservados /
NuFACE® and FIX MicroWand® are registered trademarks of the Carol Cole Company
Medical Device Safety Service GmbH